330mm

278mm

353mm

iTeraCare

Off

High Speed Low Heat

Low Speed Low Heat

High Speed High Heat

Weight loss and slimming

Anti - aging

Reduces blood pressure

Anti - inflammation

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iTera - Premium

5TH GENERATION

- INCREASES ALERTNESS LEVEL
- FUELS THE BRAIN
- REDUCES SICKNESS
- STRENGTHENS THE IMMUNE SYSTEM
- REDUCES TOXINS IN THE BODY
- HYDRATES THE BODY ULTRA FAST
- REINFORCES HEALTHY WEIGHT LOSS
- IMPROVES COMPLEXION AND SKIN RADIANCE
- NEUTRALIZES FREE RADICALS
- INCREASES CELLULAR OXYGENATION

BENEFITS OF DRINKING

iTeraCare™

WATER

BEST WHEN

Warm

charge your water! and enjoy these amazing benefits!

BDO
GLOBAL

Steps of using iTeraCare™

Warm up the body. Blow step 1 to 10, each part for about 2-3mins

 <p>Drink 2 glasses of warm water</p>	 <p>1. Palms</p>	 <p>2. Soles</p>
 <p>3. Behind the Knees</p>	 <p>4. Groin</p>	 <p>5. Armpit</p>
 <p>6. Lymph</p>	 <p>7. Upper Back</p>	 <p>8. Head</p>
 <p>9. Butt</p>	 <p>10. Blow Downwards</p>	 <p>Drink 2 more glasses of warm water</p>

iTeraCare device can also be blown on Head, Face, Body, Limbs and Pain areas.

*High Blood Pressure user don't blow no. 8


iTeraCare™

Enabling you to put the power
of healing into your own
hands naturally

TERAHERTZ4U

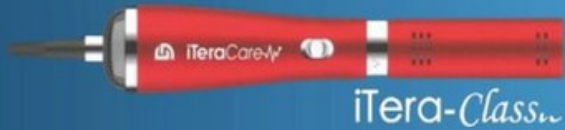iTera-Class

Prife
International

BUSINESS PLAN

iTeraCare

RECOMMENDED PROTOCOLS

iTera - *Bio*

iTera - *Bio Lite*

· Drink 1-2 glasses of warm water (preferably iTeraCare charged water) before using iTera-Bios.

· Avoid direct fan or air-conditioner blowing while using iTera-Bios.

· Please remove all metallic objects from your body during use (ie. Watch, bracelet, necklace, etc.).

· Avoid metal contact with any part of your body, hands, or legs during use.

· Ensure no usage of mobile phones or other electrical devices during use.

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iTera - *Bio*

**Clears blockages,
enhances circulation from inside out.**

Experience the resonance effects of
Terahertz waves, Bio-magnetic energy
and human cells all at the same frequency.

Dispels cold
and dampness

Promotes
circulation

Cells
activation

Diathermy
fat burner

Promotes
metabolism

Body shaping
& anti-aging

Purifies
the body

Relieves
inflammation

· Avoid touching your knees or ankles during use.

· For first-time users it is advisable to adjust the strength to below 5. You may adjust accordingly after a few sessions, depending on your health condition and body reaction.

· Drink another 1-2 glasses of warm water (preferably iTeraCare charged water) after iTera-Bio use.

· Rest for at least 1 hour before taking a shower. And avoid direct fan or air-conditioner blowing immediately after use.

· One iTera-Bio session is automatically set to 30 minutes by default. You may end the session anytime earlier before that. Let your body adjust accordingly based on your health condition. Starting with strength 1-2 and ending in 10-20 minutes is good practice for first time users.

· Start with 1 session per day for the first week.

· On the second week you may do 2 sessions per day. Ensure there is a rest period of at least 4-6 hours between each therapy. Limit to maximum of 2 sessions per day.

· It is good practice to do a session with a higher strength in the morning (ie. 4-5) and reduce to a lower strength (2-3) in the evening. This is to tone down your body before sleep.

· Do not use this appliance near baths, showers, washbasins, or other containers containing water.

· Do not touch the working parts of the appliance with metallic objects.

· During operation, both feet should fully contact with the energy board. Single point contact with the energy board is prohibited to avoid burns.

· Do not place the appliance on an uneven place such as on a slope.

THE TERAHERTZ, PEMF, & FAR INFRARED

Generated by *iTera-Bio* & *iTera-Bio Lite*

Penetrate 2-3 inches into your feet which is sufficient to **energize** your feet **acupuncture points** and your **blood circulating** to every cell in your body about every **minute**.

Prife
International E.U.

DISCLAIMER:
iTeraBio & iTeraBio Lite are not medical devices. The manufacturer does not assert its ability to treat or cure any medical conditions. Consult your medical professional for advice on any health concerns.

THE MYSTERIOUS MELODY BETWEEN HEAVEN AND EARTH

Terahertz Waves

Did you know? The frequency of terahertz waves is one trillion times per second, located in the middle of the electromagnetic spectrum between microwaves and infrared.

The non-interfering terahertz waves from nature

Every night, we bathe in non-interfering terahertz waves. These waves come from the moon, stars, and even the dark corners of the universe. Even though we can't directly perceive them, they have been around us for billions of years.


Terahertz waves and our lives

Terahertz waves have a wide range of applications, From enhancing wireless communication speed, to new medical imaging technology, to material analysis and exploring the mysteries of the universe.

These possibilities will profoundly impact our lives in the future.

LET'S STEP INTO THE WORLD OF TERAHERTZ WAVES.
WE ARE EXPLORING ENDLESS POSSIBILITIES.

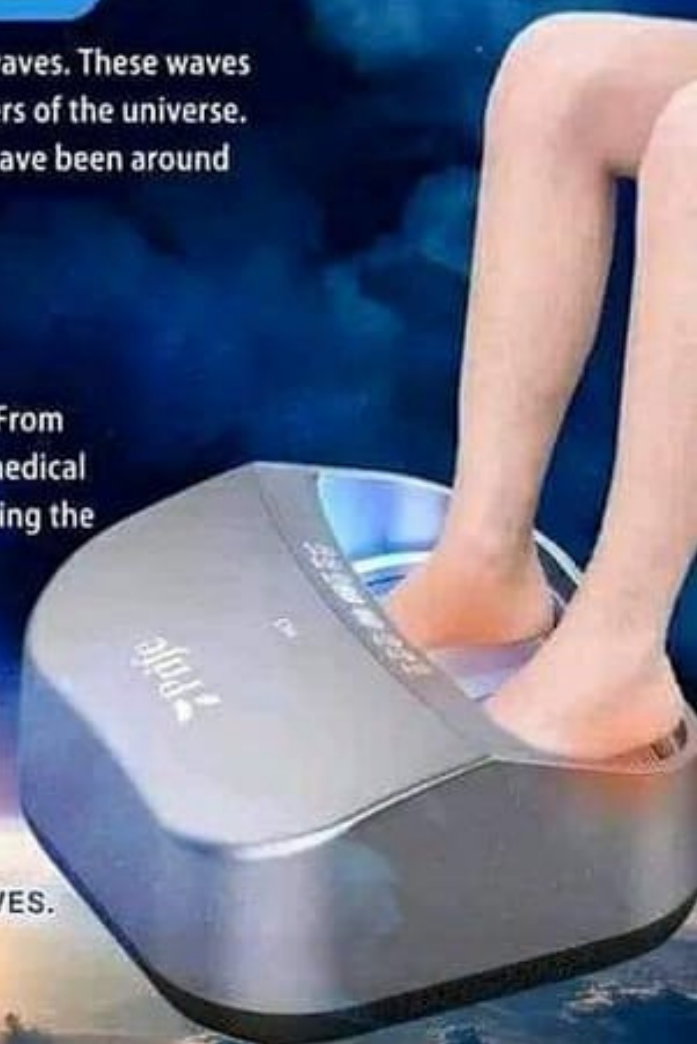Double Energy for Better Effect

Energy 1

Releases 8-1000 μ m precise Terahertz waves compatible with the human body to replenish cellular energy.

Deep penetration into the subcutaneous up to 4-5cm, which is safe and harmless.

Energy 2

Releases 1MHz Bio-RF Magnetic Energy
Improves circulation,
stimulates cellular immunity.

Improves circulation
Stimulates cellular immunity

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PRIFE IONShield

Wear inside or outside clothing

60 million per cm³ negative ions, a highly concentrated filtration

Continuous absorption into our body which is good for health, increases metabolism, improves digestive system, blood circulation and more.

O₂


NEGATIVE IONS REPEL
ACQUINAMENTS & POLLUTION
FROM AFFECTING OUR HEALTH

Prife

IONShield
Ionizer Pendant

Clear
White

IONShield

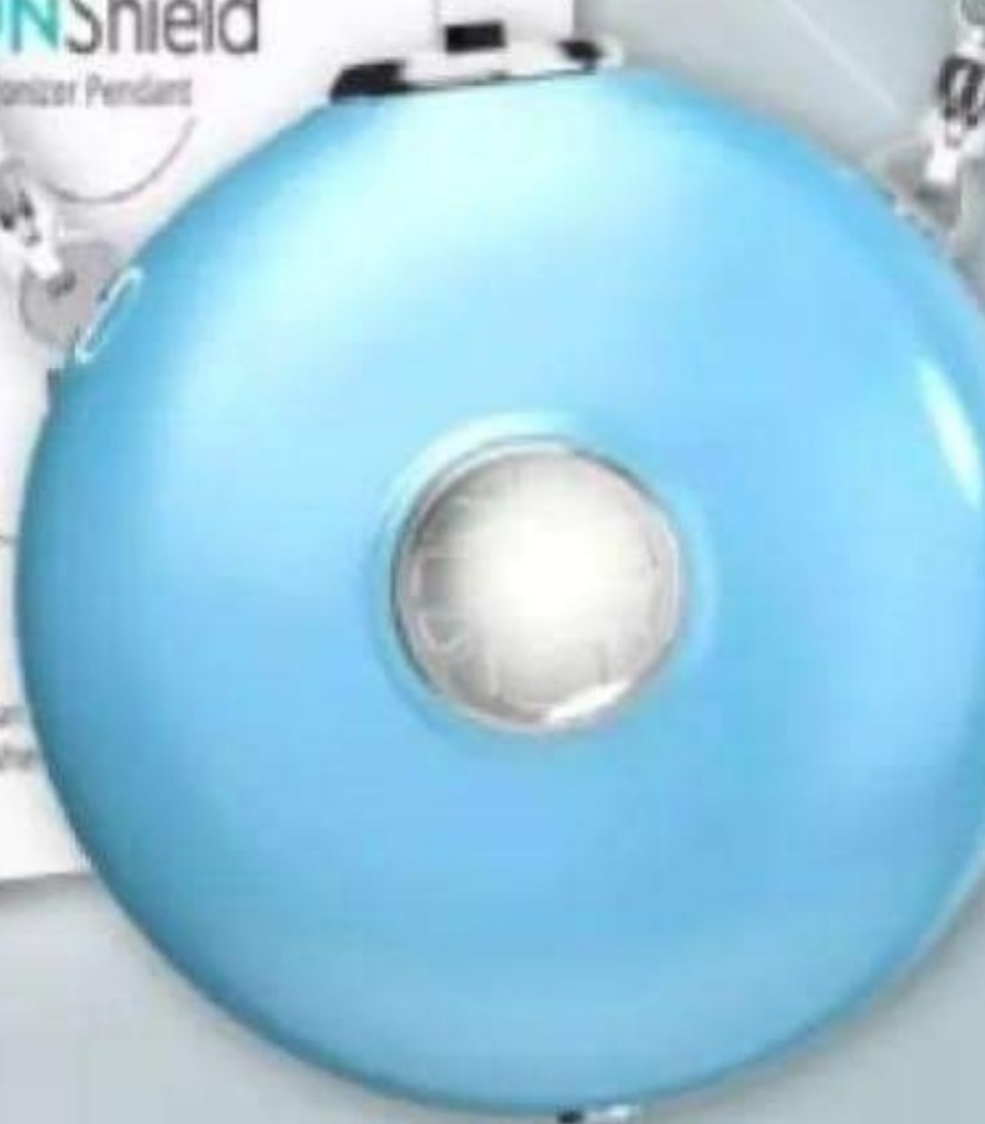“Future medicine will be the medicine of frequencies.”

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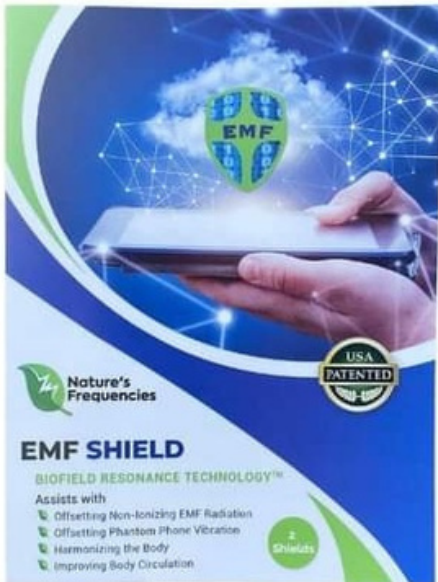
“If you want to find the secrets of the Universe, think in terms of energy, frequency and vibration.”

— Nikola Tesla

Microwave Cellphone Effects

Absorption in the Brain According to Age


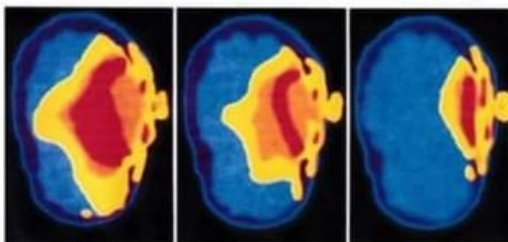
Image courtesy of Dr. Om Gandhi, University of Utah, 1996, IEE Publication

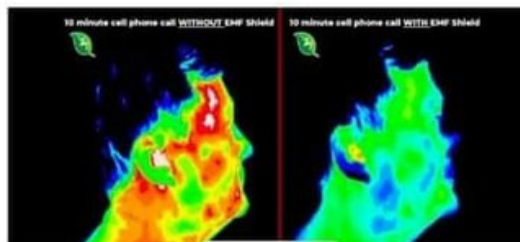Skann meg!

“Future medicine will be the medicine of frequencies.”

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"If you want to find the secrets of the Universe, think in terms of energy, frequency and vibration." – Nikola Tesla

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